

MORTECHTM
MANUFACTURING CO., INC.
ISO 9001:2015 CERTIFIED COMPANY

Family
Owned &
Operated

Celebrating
over
35
years
of service



Anatomical Pathology
Equipment Manufacturer

Elevating Pathology Workstation with Down-draft and Back-draft

GL100 Pathology Workstation

Model GL100 Elevating Pathology Grossing Station is our premiere unit. It comes standard with all required features necessary for any grossing procedure. We carry a complete listing of optional accessories to fulfill any facility requirement.

Standard Design Features

- All stainless steel construction of the entire work area and all exterior closure panels.
- Height adjustment from 32 in./81 cm. to 44 in./111 cm. via electro-hydraulic mechanism
- Exhaust duct for connection to building ventilation system with adjustable bellows
- Unit features both down-draft and back-draft ventilation
- Foot operated faucet control for hot and cold water
- Spray hose assembly with hand control
- Heavy-duty commercial waste disposer, 1/2 hp and switch
- Dissecting area rinse to provide a constant flow of water below the work area
- Top mounted LED light fixture
- Stainless C-fold paper towel holder
- Magnetic instrument holder
- Polyethylene dissecting board (Photo Blue)
- GFCI duplex receptacle with waterproof cover



GL100

Comes standard with all required features for any grossing procedure



Mortech Manufacturing Co., Inc.
411 North Aerojet Avenue
Azusa, CA 91702
Toll Free: 800-410-0100
www.mortechmfg.com

Elevating Pathology Workstation with Back-Draft Ventilation

GL105



GL105 Pathology Workstation

Model GL105 Elevating Pathology Gross Station is our elite unit. All required features for any grossing procedure comes standard. We carry a complete listing of optional accessories to fulfill any facilities requirements.

Standard Design Features

- All stainless steel construction of the entire work area and all exterior closure panels
- Height adjustment from 32 in./ 81 cm. to 44 in./111 cm. via electro-hydraulic mechanism
- Exhaust duct for connection to building ventilation system with adjustable bellows
- Large integral sink with mixing faucet
- Foot operated faucet control for hot and cold water
- Spray hose assembly with hand control
- Heavy-duty commercial waste disposer, 1/2 hp and switch
- Dissecting area rinse to provide a constant flow of water below the work area
- Top mounted LED light fixtures
- Formalin container with spigot
- Stainless C-fold paper towel holder
- Magnetic instrument holder
- Polyethylene dissecting board (Photo Blue)
- GFCI duplex receptacle w/ waterproof cover
- Length available: 48 in., 60 in., 72 in.



Countertop Pathology Grossing Station



MORTECH Grossing stations are standard type 304 stainless steel with #4 finish. For an added charge type 316 is available.

GL110 Countertop Pathology Workstation

An economical counter-top pathology workstation with all the necessary fixtures and accessories to perform routine grossing procedures.

MODEL GL 115 Base Leg Assembly


Standard Design Features

- Available in two sizes: 48 in./121 cm. long and 60 in./152 cm. long
- All stainless steel construction of the entire work area and all exterior closure panels.
- Integral sink with mixing faucet
- Exhaust duct for connection to building ventilation system
- Top mounted LED light fixtures
- GFCI duplex receptacle with waterproof cover
- (2) incandescent lights with 3x magnifier mounted on flexible arms
- Formalin container with spigot
- Stainless C-fold paper towel holder
- Magnetic instrument holder
- Polyethylene dissecting board (Photo Blue)



Mortech Manufacturing Co., Inc.
411 North Aerojet Avenue
Azusa, CA 91702
Toll Free: 800-410-0100
www.mortechmfg.com

Pathology Grossing Station

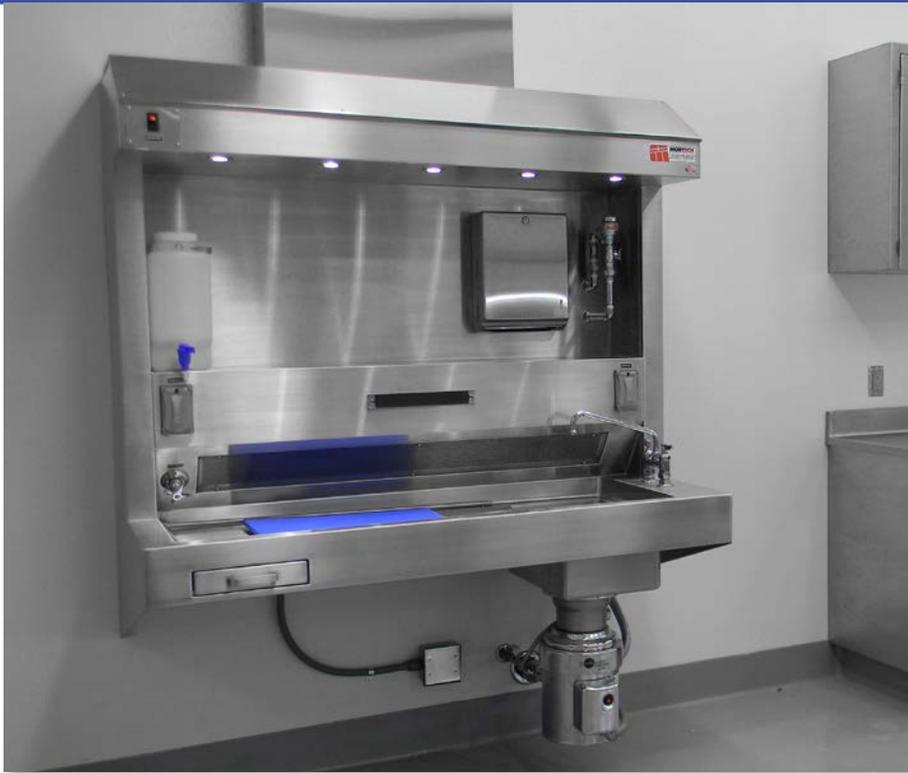

**Wall Mounted
GL115**

**Model GL125
Base Assembly
with Cabinet**



GL 115 Pathology Workstation

The Model GL115 has the same useful features as the Model GL110 with the added bonus of a larger sink basin. The larger sink basin allows the unit to be equipped with an optional disposal unit.

Standard Design Features

- Available in two sizes: 48 in./121 cm. long and 60 in./152 cm. long
- All stainless steel construction of the entire work area and all exterior closure panels.
- Stainless steel base frame with adjustable feet
- Integral sink with mixing faucet
- Exhaust duct for connection to building ventilation system with adjustable bellows
- Top mounted LED light fixture
- Incandescent light with 3x magnifier mounted on flexible arm
- Formalin container with spigot
- C-fold paper towel holder, stainless steel
- Magnetic instrument holder
- Polyethylene dissecting board (Photo Blue)



Pathology Grossing Stations

Options & Accessories

Options based on available space

When planning an order, we will assist you in selecting options to complete your system. Feel free to contact your sales associate.



OA-103 - HANDS FREE CONTROLS - HOT/COLD WATER Individual control for easy on/off activation of hot/cold water. Hot and Cold water sensors always sold together.



SENSOR - HANDS FREE CONTROLS



OA-105 - HANDS FREE CONTROLS - DISPOSAL Individual control for easy on/off activation of disposer.



OA-106 - HANDS FREE CONTROLS - END RINSE Individual control for easy on/off activation of end rinse.



OA-107 - HANDS FREE CONTROLS Individual control for easy on/off activation of table rinse, mixing faucet, waste disposal, and hot/cold water.



OA-142 - ARTICULATING VIDEO CAMERA MOUNT Articulated arm with large locking knob with variable friction on an adjustable slide track. Supplied with a quick release camera plate and secondary safety lock.



OA-110 - ARTICULATING CAMERA MOUNT Articulated arm with large locking knob with variable friction on an adjustable slide track. Supplied with a quick release camera plate and secondary safety lock.



OA-109 - MICROPHONE ON FLEX ARM Easy-to-adjust, small-diameter, alternating gooseneck provides a steep low frequency attenuation to improve sound pickup without affecting voice quality.



OA-134 - VOICE ACTIVATED RECORDER Easy to use voice recorder. Equipped with a front facing LCD screen, both convenient and compact.



OA-134-FP - FOOT PEDAL - VOICE ACTIVATED RECORDER Foot pedal allows for hands-free control of the recording and playback process.



OA-135 - KEYBOARD MOUNT ONLY
Allows full articulation of the keyboard, extends up to 11 in. Available as wall mount or 19 in. or 36 in. track.



OA-136 - MONITOR AND KEYBOARD MOUNT Allows full articulation of the flat panel monitor and keyboard, folds up over the mount and extends up to 11 in. Available as wall mount or 19 in. or 36 in. track.



OA-137 - MONITOR AND KEYBOARD MOUNT - SEPARATE Allows full individual articulation of the flat panel monitor and keyboard, folds up over the mount and extends up to 11 in. Available as wall mount or 19 in. or 36 in. track.



OA-138 - DUAL MONITOR AND KEYBOARD MOUNT Holds (2) monitors horizontally or vertically on a fully articulating arm with mouse and keyboard supports. Folds up over the mount and extends up to 11 in. Available as wall mount or 19 in. or 36 in. track.



OA-164 - SINGLE MONITOR MOUNT Allows full articulation of the flat panel monitor, folds up right over the mount and extends up to 11 in. Available as wall mount or 19 in. or 36 in. track.



OA-165 - DUAL MONITOR MOUNT Holds (2) monitors horizontally or vertically on a fully articulating arm. Folds up over the mount and extends up to 11 in. Available as wall mount or 19 in. or 36 in. track.



OA-166 - SINGLE TABLET MOUNT
Spring loaded cradle expands and contracts to fit most tablets. Securely holds with exceptional range of motion.



OA-167 - LAPTOP MOUNT
Allows full articulation of Laptop platform, folds up over the mount and extends up to 11 in. Available as wall mount or 19 in. or 36 in. track.

Build Your Own Mount

Mortech can offer numerous configuration options, solutions and combinations:

- Over 2,000 modular arm designs that can be adapted to your needs.
- Custom bracket and attachment options to fit your specific device.
- Fully custom mounts to meet your unique challenge.

When configuring your custom setup, please consider these questions:

- *Mounting to what surface?*
- *Mounting a flat panel monitor?*
- *Mounting a keyboard?*
- *Mount type?*
- *What arm segments?*
- *Arm segment length?*
- *Separate arms?*
- *Keyboard tray size?*
- *Keyboard tray type?*




OA-139 - STATIONARY MONITOR MOUNT ON BACKSPASH Flat panel monitor attaches directly to backplash of grossing station.



OA-145 - CPU WALL MOUNT
Designed to fit CPU's from 4.6 to 8 in. wide. Units are constructed from heavy-duty aluminum, with neoprene padding to protect your CPU as it is clamped securely into the mount.


Pathology Grossing Stations

Options & Accessories

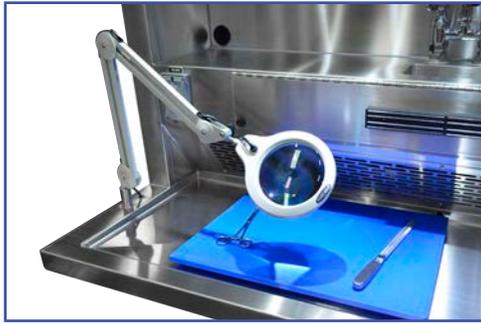

OA-124 - LED MAGNIFIER LIGHT

Flexible spring-balanced arm, handle in shade permits easy positioning.

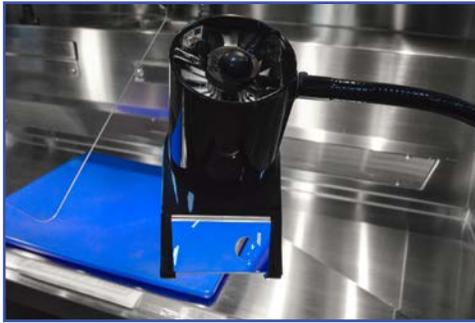

OA-123 - FLEX ARM MAGNIFIER LIGHT

Flexible spring-balanced arm, handle in shade permits easy positioning. Retractable magnifying extension for closer examination.



OA-156 - FLEX ARM LED LIGHT

Adjustable to all angles using a 30 in. flexible arm and powerful LED bulb.



SLSE50-DK - SURGICAL EXAM LIGHT

Delivers best in class intensity and color. Dual lens optical system collects and concentrates light into a precise beam.



M184042 - SLMAG - MAGNIFYING EXAM LIGHT

World class optics, performance and precision movement. The perfect combination of light quality, positionability and large viewing area.




WAVE - LED MAGNIFIER LIGHT

Large viewing area and LED energy saving technology with a 3.5 diopter magnifier.



OA-130-110 - SINGLE X-RAY ILLUMINATOR

Available in 110 or 230v. "Instant on" lighting, maximum quality and flexibility for radiological film interpretation. Ideal for use in hospital and medical high-volume environments.

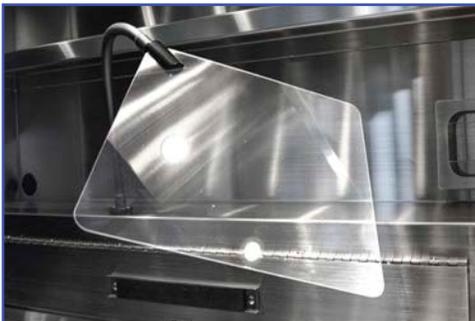

OA-141 - SAFETY SPLASH GUARD

Mounted per customer specification. Articulating arm allows protection for both specimen and user from cross contamination.



OA-146 - FULL PLEXIGLASS ENCLOSURE

Covers and protects the working area for a clean and sanitary area. Constructed of 0.5 in./1 cm. plexiglass, contains fumes and reduces contaminants.



OA-155 - DRY ERASE BOARD

Provides easy dry-erase writing for meetings, brainstorming and discussions. Strong fiberboard and frame resist warping and bending.



OA-144 - MAGNETIC INSTRUMENT HOLDER

Conveniently store and organize your instruments.




OA-125-24, OA-504-48 - REMOVABLE MEASURING RULER

Available in 24 in. and 48 in.. Provides a convenient measuring device, easily accessed and mounted in plain sight. Also available in 48 in. Provides a convenient measuring device, easily accessed and mounted in plain sight.



OA-147-1, OA-147-2, OA-147-3 GLOVE BOX HOLDER Stainless steel support to hold latex and nitrile glove boxes. Holders come in a side or top load model. Available in 1, 2, & 3 capacity.

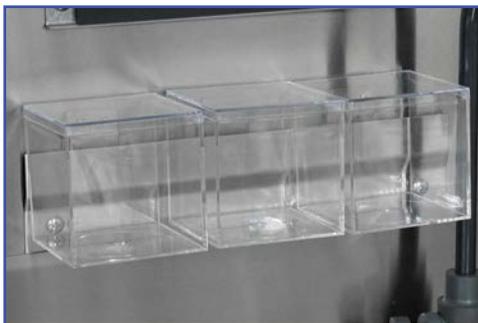

OA-125 - CASSETTE HOLDERS - SET OF 3 Conveniently holds cassettes within arms reach and removable boxes mount to rail on front of unit, Set of (3) boxes provided



OA-132 - STAINLESS STEEL STORAGE CONTAINERS Set in a No. 4 finish. Removable containers come (1) large and (1) small each, per order.



OA-158 - STAINLESS C-FOLD PAPER TOWEL HOLDER Surface mounted towel dispenser. Fabricated of type 304, 22 gauge stainless.



OA-159 - STAINLESS C-FOLD PAPER TOWEL HOLDER - OPEN DESIGN Surface mounted towel dispenser with open design.



OA-148-1 - ADJUSTABLE STAINLESS STEEL SHELF - ONE All stainless steel construction adds storage space for accessories. Adjustable shelf is 12 in. W x 7.5 in. L.



OA-148-2 - ADJUSTABLE STAINLESS STEEL SHELF - TWO All stainless steel construction adds storage space for accessories. Adjustable shelf is 12 in. W x 7.5 in. L.



OA-149 - STATIONARY STAINLESS STEEL SHELF - ONE All stainless steel construction adds storage space for accessories. Stationary shelf is 12 in. W x 7.5 in. L.



OA-133 - FORM HOLDER Sturdy construction and can be mounted to store documents away from any fluids and risk of damage, (2) provided standard.



OA-127 - UTILITY DRAWER Constructed of all stainless steel, this drawer is ideal for additional storage of utensils, office supplies or personal items



OA-128 - WRITING PLATFORM Sturdy stainless steel construction provides a reliable smooth surface for writing and conveniently hides away.



OA-129 - COMBO UTILITY DRAWER / WRITING PLATFORM Conveniently combines both accessories in one, a writing platform with storage drawer underneath.

Pathology Grossing Stations

Options & Accessories



LW418 - ELECTRONIC SCALE - PLATFORM
6.5 x 5.5 in. Platform, 2000 gm.
 Compact, lightweight, portable and perfectly suited for laboratory and educational purposes. 2200 x 1 gm.



LW419 - ELECTRONIC SCALE - PLATFORM
6.6 x 6 in. Platform, 2100 gm.
 Compact, lightweight, portable and perfectly suited for laboratory and educational purposes. 2200 x .1 gm.



LW425 - ELECTRONIC SCALE - PLATFORM
6 in. Platform, 2000 gm.
 Compact, lightweight, portable and perfectly suited for laboratory and educational purposes. 2200 x .1 gm.



OA-154 - SHARPS DISPOSAL CONTAINER
 Solution for disposal of needles, syringes, and lancets, as well as other types of sharps and small quantities of medical waste.



OA-153 - SPECIMEN RINSE BASKET/HOLDER FOR SINK Single mesh fine with a plastic handle includes stainless holder that fits over sink.



OA-152 - HANDS FREE SOAP DISPENSERS
 High-efficiency pump mechanism provides quick and precise dispensing. Easy-to-use button control and the wide opening makes refills fast, easy and mess-free. (2) models available: Deck and Wall Mount.



OA-126 - VENTILATED TRASH CONTAINER
 Stainless steel construction, removable top lid, holds 13 gallon trash liner, 1.5 in. duct attaches to station exhaust system for fume removal, Container size: 13 in. x 13 in. x 29 in.



OA-140 - FLOOR MOUNT ANCHORS
 Securely holds your grossing station in place with sturdy gauge floor mount brackets. Engineered seismic calculations available at an additional cost.



OA-161 - STAINLESS STEEL TISSUE MARKING DYE TRAY Conveniently stores tissue dye bottles within reach and easily accessible.



OA-108 - SELF CONTAINED VENTILATION ASSEMBLY Consisting of 3 blowers with filters to clean and safely recycle air within your work area.



OA-143 - REPLACEMENT EXHAUST BELLOWS
 Adjustable bellows replacement for exhaust duct connection to building ventilation system.



GL113 - REPLACEMENT FILTERS - MODEL GL113 Activated filter for odorous and formaldehyde fumes.



OA-100-DM - CO2/NATURAL GAS VALVES
All brass with chrome plated finish. Deck mounted. Visually marked for type of function required.



OA-101-DM - COMPRESSED AIR VALVE
All brass with chrome plated finish. Deck mounted. Visually marked for type of function required.



OA-102-DM - VACUUM VALVES
All brass with chrome plated finish. Deck mounted. Visually marked for type of function required.



OA-100-WM - CO2/NATURAL GAS VALVES
All brass with chrome plated finish. Wall mounted. Visually marked for type of function required.



OA-101-WM - COMPRESSED AIR VALVE
All brass with chrome plated finish. Wall mounted. Visually marked for type of function required.



OA-102-WM - VACUUM VALVES
All brass with chrome plated finish. Wall mounted. Visually marked for type of function required.



OA-111 - EMERGENCY EYEWASH / HOSE ASSEMBLY
Hand controlled spray to provide a gentle flow of water to rinse eyes of foreign debris.




OA-112 - EYEWASH / FAUCET ASSEMBLY
Deck mounted eyewash for hands free eye wash or can be removed for use as a drench hose for eye, face, or body use coverage.



OA-113 - EYEWASH - FLIP DOWN STYLE
Wall mounted polished chrome-plated brass single action pull down activated eye/face wash.



OA-150 - SPRAY HOSE ASSEMBLY
Easy grip assembly with flexible hose, conveniently placed for easy spray cleaning.

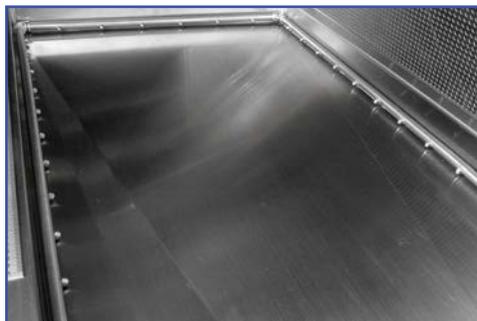

OA-115 - THREE SIDED RINSE ASSEMBLY
Washes debris towards sink basin controlled by cold water valve or hands free control.



OA-116 - END RINSE ASSEMBLY
Washes debris towards sink basin controlled by cold water valve or hands free control with back flow prevention.

Pathology Grossing Stations

Options & Accessories




LC-50 - 1/2 HP HEAVY DUTY COMMERCIAL DISPOSER

Standard on most sinks. Variable power options. Available in 1/2, 1 or 2 HP.



SS-100, SS-200 - HEAVY DUTY COMMERCIAL DISPOSER

Variable power options. Available in 1 or 2 HP.

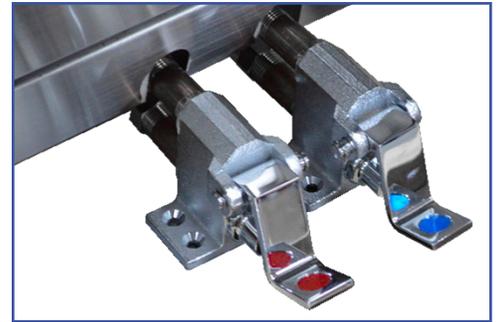

OA-151 - HOT AND COLD FOOT VALVE

Sturdily constructed hot/cold foot pedal helps temper water and keep your hands free while working.



OA-117 - FLAMMABLE STORAGE CABINET

NOT FOR ALL UNITS - Protect workers, reduce fire risks, and improve productivity by storing hazardous liquids

Options based on available space

When planning an order, we will assist you in selecting options to complete your system. Feel free to contact your sales associate.



OA-118 - FORMALIN DISPENSING AND COLLECTION

2.5 gal. container and receptacle hold reserve formalin and collects waste. Foot controlled motor dispenses formalin at a steady flow.



OA-157 - FORMALIN CARBOY W/ SPIGOT


Space saving rectangular shape and chemical resistant holds 9 liters.

Multiple Formalin Options Available:

- OA-119 - FORMALIN DISPENSING ONLY - FROM BUILDING
- OA-120 - FORMALIN DISPENSING ONLY - SELF CONTAINED
- OA-121 - FORMALIN COLLECTION ONLY - TO BUILDING
- OA-122 - FORMALIN COLLECTION ONLY - SELF CONTAINED
- OA-162 - FORMALIN DISPENSING & COLLECTION - CUBE STYLE
- OA-163 - FORMALIN DISPENSING ONLY SYSTEM - CUBE STYLE

Pictured accessories may differ depending on your unit and other optional accessories. Please review the model's standard features to ensure accuracy of specifications.



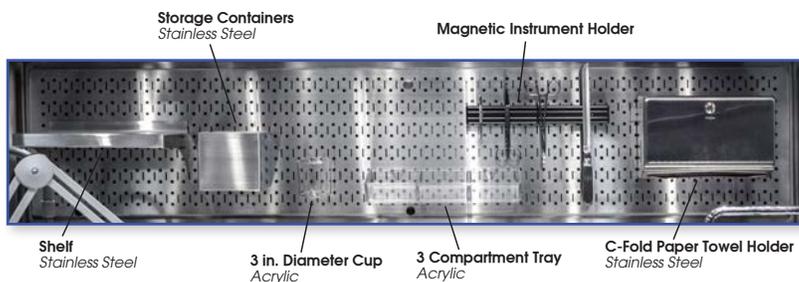

Grossing Station Personalization


Options based on available space

When planning an order, we will assist you in selecting options to complete your system. Feel free to contact your sales associate.

OA-170 - Pegboard Backplash for GL100 or GL105



OA-170-1 - MAGNETIC INSTRUMENT HOLDER Conveniently store and organize your most commonly used instruments within arms reach.

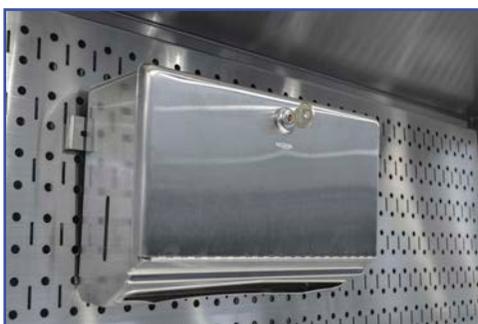

OA-170-2 - C-FOLD PAPER TOWEL HOLDER - ENCLOSED Stainless steel surface mounted towel dispenser.




OA-170-3 - C-FOLD PAPER TOWEL HOLDER - OPEN DESIGN Stainless Steel surface mounted towel dispenser with open design.



OA-170-4 & 170-5 - GLOVE BOX HOLDER Stainless Steel support to hold latex and nitrile glove box. Available in single & double box holders.



OA-170-6 - ADJUSTABLE STAINLESS STEEL SHELF All stainless steel construction adds storage space for accessories. 12 in. W x 7.5 in. L.

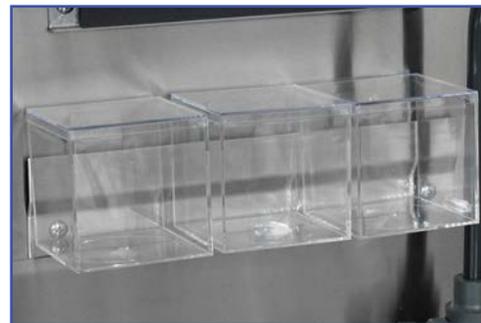


OA-170-7 - CASSETTE HOLDER - ACRYLIC Conveniently holds cassettes within arms reach. Removable boxes mount to pegboard on front of unit. Set of (3) boxes provided.



OA-170-8 - CASSETTE HOLDER - STAINLESS STEEL - SMALL Stainless steel storage container conveniently holds cassettes within arms reach. Removable boxes mount to pegboard on front of unit. 4 in. L x 3 3/4 in. H x 3 5/8 in. D.



OA-170-9 - CASSETTE HOLDER - STAINLESS STEEL - LARGE Stainless steel storage container conveniently holds cassettes within arms reach. Removable boxes mount to pegboard on front of unit. 5 in. W x 4 3/4 in. H x 4 5/8 in. D.

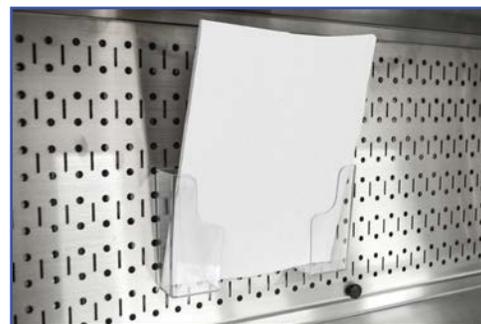


OA-170-10 - ADJUSTABLE FORMS HOLDER - ACRYLIC For the storing of documents away from any fluids and risk of damage. Adjustable, to hold any media size. 2 in. W x 5 in. H x 1.5 in. D.

PATHOLOGY GROSSING STATIONS


OPTIONS & ACCESSORIES




OA-170-11 - CLEAR CUP - ACRYLIC
Cup display for pegboards.
3 in. W x 3 in. H.

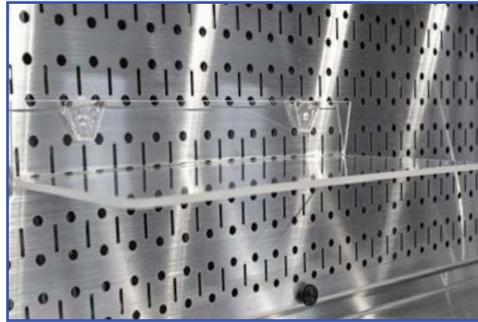

OA-170-12 - CLEAR SHELF - ACRYLIC
Clear Shelf and attached molded peg hooks for easy display.
13.5 in. W x 6 in. D.

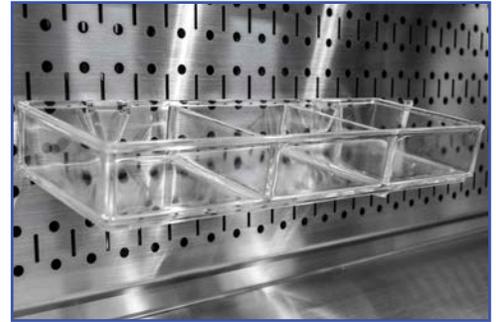

OA-170-13 - THREE COMPARTMENT CLEAR TRAY - ACRYLIC Molded tray with dividers, and attached molded peg hooks.
12.75 in. W x 1.75 in. H x 4.5 in. D.

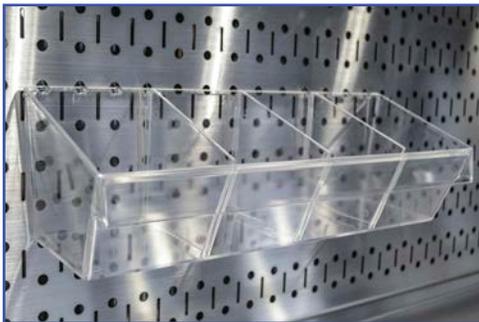

OA-170-14 - THREE COMPARTMENT CLEAR TRAY - ACRYLIC Designed for displaying multiple products of varying sizes at one time. Add included 4 in. dividers (qty. 3) to create more compartments.
13 in. W x 2.5 in. H x 4 in. D.



OA-170-15 - TWELVE COMPARTMENT CLEAR TRAY - ACRYLIC Multifunctional tray for pegboard or counter-top display. 12 separate compartments on 3 levels.
11.25 in. W x 9 in. H x 9 in. D.




OA-170-16 - 4 .IN SQUARE CLEAR BOX - ACRYLIC For the display and storage of accessories and small objects.
4 in. L x 4 in. H.



OA-170-17 - 6 .IN SQUARE CLEAR BOX - ACRYLIC For the display and storage of accessories and small objects.
6 in. L x 6 in. H.



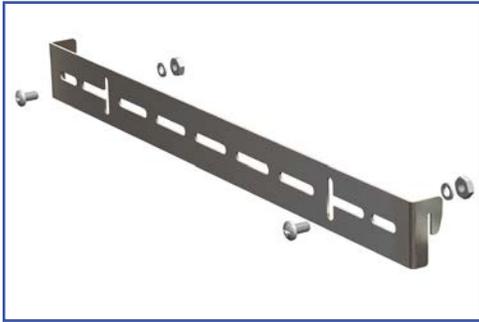

OA-170-18 - UNIVERSAL PEGBOARD MOUNTING BRACKET Provides the ability to mount any item to pegboard.



OA-154 - SHARPS DISPOSAL CONTAINER
 Pegboard mounted container for disposal of needles, syringes, and lancets, as well as other types of sharps and small quantities of medical waste.



OA-139 - MONITOR PEGBOARD MOUNT
 Attach flat panel monitor directly to pegboard

Pegboard backslash removal procedure for cleaning purposes:




Fastening knobs (2) located at lower right & left




Unscrew knobs (2)



Using handles (2) located on right & left, lift & pull to remove backslash

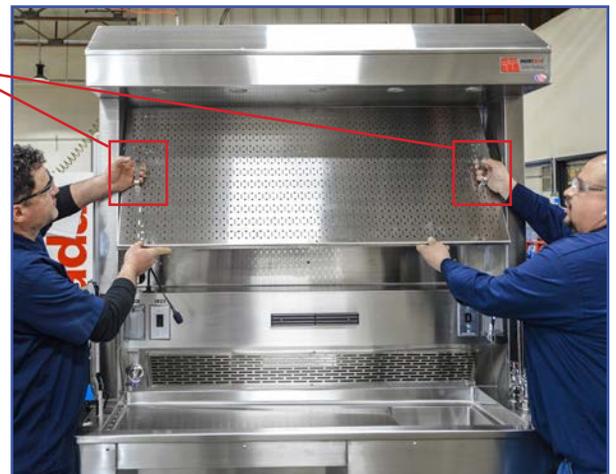

We recommend working with a partner when installing or removing backslash

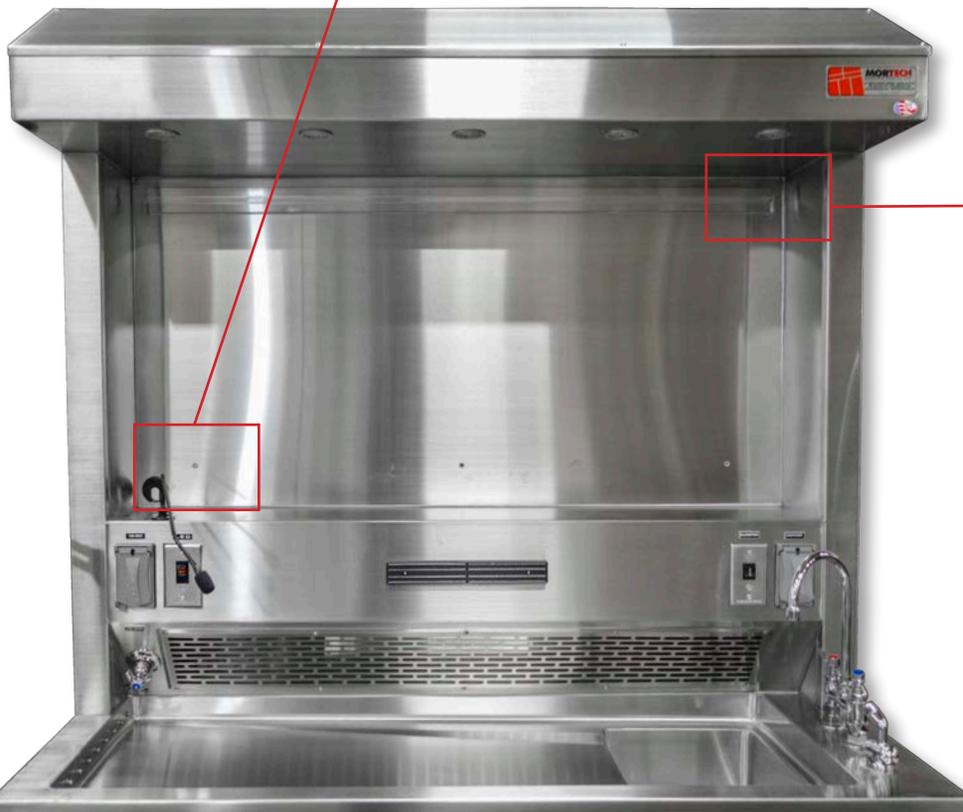

Mounting bracket located behind backslash secures pegboard to station

Grossing Station Personalization

Options based on available space

When planning an order, we will assist you in selecting options to complete your system. Feel free to contact your sales associate.

OA-170 - Pegboard Backslash for GL100 or GL105

Mobile Cabinets & Carts

Mobile Utility Cabinet

The stainless steel mobile cabinet offers a wide range of storage options. All exposed hardware is brushed stainless steel. Hinges are thick stainless steel laboratory/hospital type, five knuckle frictionless and are long with fast pin and rounded ends. All doors are provided with ball catch assembly for closing. Pulls are available in two styles, a satin chrome plated extended handle or a bar pull with a satin finish. Both are mounted to the inside head. When doors are supplied with a lock, the lock will be a removable core for enhanced security and quick rekeying. In cases of double doors, the left door will be fitted with an astragal.



Stainless Steel Utility Carts

The Stainless Steel Utility Cart has many valuable functions in either a lab or morgue setting environments. It has the ability to transport all types of laboratory ware and can also be used as an instrument/dissecting tray.

Standard Design Features

- All welded construction (except casters)
- Stainless with premium #4 polished shelves in durable 16 gauge
- Bolt on casters, 2 swivel & 2 rigid, for superior tracking
- Tubular handle with smooth radius bend-standard
- 1.5 in./4 cm. shelf lip for retention
- Top shelf height: 35 in./89 cm.
- Clearance between shelves: 12 in./30 cm.

*All measurements
in inches/centimeters +-*



Mortech Manufacturing Co., Inc.
411 North Aerojet Avenue
Azusa, CA 91702
Toll Free: 800-410-0100
www.mortechmfg.com



These ventilated storage cabinets can be ducted to the exhaust system of your facility. (10) Drawers per unit provides more than enough space to store your items. Each drawer is equipped with a lock and key.

- (5) Lockable Drawers with maximum capacity of 250 lbs.
- (4) Vented - Morgue Storage Cabinets.
- 75 to 150 CFM depending on size of the cabinet.
- 24 in./60 cm. Depth. 37 in. Standard base.

Ventilated Base Cabinets



Customized Stainless Steel Casework & Table Assemblies

Specimen Photo Table

- Column Custom: black anodized aluminum extrusion.
- Carriage: Friction drive, dual drive rollers, dual guide rollers, large drive knob.
- Counterbalance: Spring Long life constant torque.
- Dual Brake Pads: Feature positive "no-shift" carriage position lock.
- Camera Support Plate: 2-1/2 in. H x 5 in. W, rubber pad.
- Camera Alignment: Up to 5-1/2 in. forward, horizontal slot for side alignment.
- Baseboard: Black powder coated steel with gray screened centering grid.
- Light output 36W each, 72W total (equivalent to approximately 200W incandescent).
- Color temperature 5200K
- Light control: Power switch for both lights, left lamp on/off switch only.
- Electrical 120/240VAC - 1.5 in./75 cm.



Complete stand with high frequency fluorescent lights and light arms, one each side. Ideal for digital imaging and specimen photography these cool full spectrum fluorescent lamps are daylight balanced at 5200K for superb color clarity and operate at 12,000Hz to eliminate flickering in video applications. Free standing light control has on/off main switch and on/off switch for left lamp to permit shadowing on raised, embossed or bas-relief subject matter. Total 72-watt output equivalent to approximately 200-watts of incandescent lighting.



Ventilated Formalin Station

VFS-SC (Self Contained)
VFS-D (Ducted to building)

The Ventilated Formalin Station is a stainless steel base for your formalin cube equipped with a blower that filters out toxic fumes. It is also equipped with sensors complete with audible and visual alarms to alert the user of hidden leaks or breaches. The station fits the standard 5 gallon formalin cube.

Audible and Visual Alarms

The station is equipped with a visual and audible alarm activated by a formalin leak. The sensors are located at the base of the station and in the drip pan resting under the spigot of your formalin container.

Ventilated Tissue Storage



Tall Vented Storage Unit Glass Door

- Five Adjustable Shelves
- Vented - Attach to Building Exhaust System
- 75 to 150 CFM depending on size of the cabinet
- Glass Hinged Doors:
SC613-36
SC613-48

Depth available: 16 in./40 cm., 18 in./45 cm.,
22 in./55 cm., 24 in./60 cm.

Standard height: 72 in./182 cm., 84 in./213 cm.

Tall Vented Storage Unit Solid Door

- Five Adjustable Shelves
- Vented - Attach to Building Exhaust System
- 75 to 150 CFM depending on size of the cabinet
- Hinged solid doors
SC612-36
SC612-48

Depth available: 16 in./40 cm., 18 in./45 cm.,
22 in./55 cm., 24 in./60 cm.

Standard height: 72 in./182 cm., 84 in./213 cm.

Mortech's Ventilated tissue storage is intended fit anywhere in your facility. It is ducted to the buildings exhaust system and vents out any fumes your storage items may have. For improved organization and storage accessibility, racks roll in and out with ease.

Ventilated Tissue Storage Cabinets

- All Stainless Steel Construction
- Vented - Attach to Building Exhaust System
- 75 to 150 CFM depending on size of the cabinet
- (4) Adjustable leveling feet
- Glass hinged solid doors

Standard Depth: 30 in./76.2 cm.

Standard height: 72 in./183 cm., 84 in./213 cm.

Mobile Specimen Shelving

- (4) Adjustable Shelves
- (4) Casters with brakes
- Open grid mesh.

Depth available: 24 in./70 cm.

Standard height 36 in./91 cm., 48 in./121 cm., 60 in./152 cm.



Dual-Draft Pathology Workstation

- Fabricated from heavy gauge 304 stainless steel with a No. 4 finish
- Helic-welded seams and joints, ground and polished to match adjoining surfaces
- Radii on all inside and outside corners for easy cleaning
- Stainless steel base frame with adjustable feet
- Faucet with hot and cold water control
- Hose rinse assembly with cold water control
- Rapid and positive drainage
- 4 in. LED recessed lighting
- GFCI duplex receptacle with waterproof cover.
- Plumbing and electrical lines factory installed and required only single point connection at the project location



Nylon Biopsy Bags

STANDARD

SMALL BAGS

Order No	Size (Description)
BBS1000-02	Small Bags/1000 per box
Width:	1.1 in./30 mm.
Length:	1.9 in./50 mm.

MEDIUM BAGS

Order No	Size (Description)
BBM200-02	Medium Bags/200 per box
Width:	1.7 in./45 mm.
Length:	2.9 in./75 mm.

LARGE BAGS

Order No	Size (Description)
BBL500-02	Large Bags/500 per box
Width:	2.9 in./75 mm.
Length:	3.7 in./95 mm.

- Reduce risk of small specimen loss
- Design allows unrestricted fluid movement around tissue
- Bags peel open for quick and easy specimen removal
- Available in Standard & Premium
- Available in three sizes

PREMIUM

SMALL BAGS

Order No	Size (Description)
BBS1000-01	Small Bags/1000 per box
Width:	1.1 in./30 mm.
Length:	1.9 in./50 mm.

MEDIUM BAGS

Order No	Size (Description)
BBM200-01	Medium Bags/200 per box
Width:	1.7 in./45 mm.
Length:	2.9 in./75 mm.



SMALL NYLON BIOPSY BAG



MEDIUM NYLON BIOPSY BAG



LARGE NYLON BIOPSY BAG

Reduce Specimen Loss with Easy Peel Design

Design allows unrestricted fluid movement around tissue. Bags peel open easily for quick and worry-free specimen removal.

